

The Bellota Sonohysterography Catheter without Cervical Cone is available in a 5 (f) size and is 350mm in length.

The catheter is to be used for sonohysterography procedure that aid in the analysis of uterine pathology. A placement sheath eliminates the need for a tenaculum. This catheter is designed for single patient use.

PRODUCT NO. UIC 5F

PACKAGING: 25 sterile catheters per box

CAUTION

Federal law restricts this device to sale by or on the order of a physician.
Do not reuse for avoiding user may be infected by the microorganism.

INDICATIONS

The Bellota Sonohysterography Catheter without Cervical Cone (Model: UIC 5F) Catheter is used to access the uterine cavity for saline infusion Sono-hysterography procedure.

CONTRAINDICATIONS

The catheter type procedures are contraindicated in pregnancies, suspected pregnancies, active pelvic infections, or recent pelvic infections, severe uterine bleeding, gynecologic malignancies. This device should not be used for intrafallopian or intratubal procedure or for in-vitro fertilization.

PROCEDURE

- Remove the catheter from the sterile pouch.
- Attach the saline source to the open luer and allow the saline to flow through the catheter

Visit our website:
www.bioteque.com
Tel 1-800-889-9008
Fax 1-408-526-9399




MANUFACTURER

Panpac Medical Corporation
6F-2, No 202, Sec. 3, Ta-Tong Road,
Shi-Chih City, Taipei Hsien, 22103
Taiwan



STERILE EO

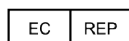

Wellkang Ltd/ The Black Church,
St. Mary's Place, Dublin 7,
D07P4AX, Ireland

CE 2460

Not made with natural rubber latex

- removing air and flushing the saline channel.
- Adjust the protective sheath so that just 1/2 to 1 cm of the catheter protrudes beyond the sheath.
- Expose the cervix using a speculum and clean it with an antiseptic solution.
- Carefully advance the 5f catheter through the cervical canal and past the internal cervical OS, so that the catheter is within the uterine cavity while the catheter is advanced.
- The 7 centimeter-marking on the catheter is helpful to place the catheter properly.
- Attach saline to open syringe and introduce through catheter slowly.
- Remove the speculum and insert the transvaginal transducer into its proper position adjacent to the catheter.
- Now activate the ultrasound system and inject saline as required by the procedure. Usually, 3cc to 20cc is sufficient saline for a proper diagnosis.
- Perform Sonohysterography in the usual way.
- Endometrial distention usually can be maintained throughout the procedure. If uterine distention is not visualized with ultrasound while saline is injected into the catheter or leaking is observed from the vagina, verify that the catheter is positioned correctly in the uterine cavity before continuing.
- Upon completion of Sonohysterography, remove the catheter carefully from uterine cavity.

STORAGE: Store at room temperature.



DO NOT REUSE



EXPIRY DATE



BATCH CODE



STERILE BY ETHYLENE OXIDE

Made in Taiwan